E07/C E07/S
E07/Ccs E07/Scs

from 04.2017

ISTRUZIONI

DIRECTIONS

MONTAGEANWEISUNGEN

INSTRUCTIONS

INSTRUCCIONES

**ISTRUZIONI
INSTRUCTIONS
ANWEISUNGEN
INSTRUCTIONS
INSTRUCCIONES**

ISTRE07

DATA: 10.03.2017 REV.4

PAGINA: 3/39

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Capitolo uno / Chapter 1 / Abschnitt 1 / Chapitre 1 / Capítulo 1 7-10

1.1 Campo di impiego
1.1 Field of use
1.1 Einsatzgebiet
1.1 Utilisations possibles
1.1 Campo de empleo

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The products shown in this instruction booklet have been designed for use in dental laboratories. Any other use must be regarded as incorrect and may be the cause of damage both to the equipment and the operator. In this case, the manufacturer does not accept any responsibility. The ideal environmental conditions for use of the equipment are:

- Indoor use
- Altitude of up to 2000 m.
- Temperature from 5 to 40°C.
- Maximum relative humidity of 80% up to 31°C in linear drop to 50% at 40°C.
- Supply voltage fluctuations +0 ±10% nominal voltage
- Transient overvoltage according to installation category **II**, 2 degrees of contamination.

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Mod. E07/C, E07/S

Suction unit for polishing motor.

This small-sized unit offers ample space for housing the polishing motor and allows the technician full ease of movement.

Lighting is provided by a sealed 18 W fixture.

A large drawer is provided for the brushes. It is quite easy to reach the suction motor and the filters, because the upper part of the unit is removable. As additional motor protection, a supplementary filter of 3654 cm² is supplied. The motor has an average life of approximately 2500 hours, and doesn't need the brushes replaced. It is necessary to check the 110V motor brushes every 700 working hours. The filter bags have a filtering area of 1500 cm².

The suction unit can be supplied in two different versions: with (mod. E07/C) or without (mod. E07/S) the polishing motor.

The power of the polishing motor is 350 W with two speeds: 2800 and 1400 r.p.m.

Mod. E07/Ccs, E07/Scs

Suction unit for polishing motor with centralized extraction.

Model identical to the previous one but set up with appropriate ducting and control microswitches for connecting to an Erio centralized extraction system.

Lighting is provided by a sealed 18 W fixture. A large drawer is provided for the brushes.

The unit can be supplied in two different versions: with (mod. E07/Ccs) or without (mod. E07/Scs) the polishing motor.

The power of the polishing motor is 350 W with two speeds: 2800 and 1400 r.p.m.

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The E07 machine must be installed by qualified personnel trained for the purpose.

Various safety rules must be observed for the safety of the machine.

- Position the polishing motor on a flat stable surface.
- Check that the supply voltage, indicated on the E07 label, corresponds to that of your mains.
- Connect the machine to a power socket with an efficient earth connection.
- Following installation, the connection plug must remain accessible to the user so that the machine can be disconnected during maintenance operations.
- The machine supply lead must be connected to the provided connectors inside the connector box (Z), Photo 1.
- Do not cover the machine's air evacuation outlets.

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-Attention! To maintain the protection level of the machine during maintenance operations, the operator must refer to the instruction manual.

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The filter maintenance operations can be carried out by the operator who uses the machine habitually provided that he/she is trained appropriately, after disconnecting the machine power supply.

Photo 1 Raise the cover of the apparatus (A) and secure it in the open position by means of the rod B. Remove the hose (C) from the intake. Remove the filter support frame (D) from the intake chamber (E), by raising it.

Open the clamp and remove the filter bag. It is recommended that the filter bag be changed every 100 effective working hours.

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2.3 Replacing the brushes of the extractor motor

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The suction motor maintenance and cleaning and/or repair operations must be carried out by qualified personnel authorized for the purpose after disconnecting the machine power supply.

Disconnect the electrical cable which powers the motor inside the box (F). By actuating the lever (G, photo 1) remove the motor box from the apparatus. Remove the cover (H, photo 2) of the motor box. Disconnect the four screws that fix the motor plate support (M). After disconnection of the electrical socket contacts (taking note of the exact position/colour of the wires) remove the plate and the motor (I). Remove the four screws that fix the motor carbon brushes (L): use a screwdriver as a lever to remove the brushes and disconnect the faston contact. Fit the new brushes, making sure that the brushes holding-cartridge fit into place, fix it with the screws and reconnect the fast-on contact. Put a little screwdriver and prize it on the upper slit of the brush cartridge to allow it to be disconnected and stuck to the motor. To complete the maintenance work, close the motor housing and carry out the disassembly instructions in the reverse order.

the faston contact. Fit the new brushes, making sure that the brushes holding-cartridge fit into place, fix it with the screws and reconnect the fast-on contact. Put a little screwdriver and prize it on the upper slit of the brush cartridge to allow it to be disconnected and stuck to the motor. To complete the maintenance work, close the motor housing and carry out the disassembly instructions in the reverse order.

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Disconnect the electrical cable which powers the motor inside the box (F). By actuating the lever (G, photo 1) remove the motor box from the apparatus. Remove the cover (H, photo 2) of the motor box. Disconnect the four screws that fix the motor plate support (M). After disconnecting the electrical socket contacts (after taking note of the exact position/colour of the wires) remove the motor (I) and replace it. Repositioning the motor, take a great of care that it adheres perfectly to the seal, if it does not there will be no suction. To complete the maintenance work, close the motor housing and carry out the disassembly instructions in the reverse order.

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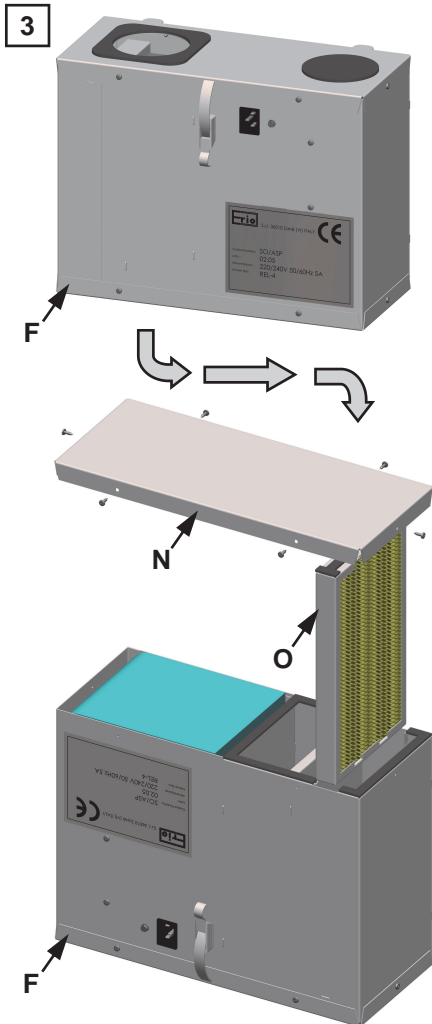
2.5 Replacing extractor motor filter cartridge

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Photo 3 Disconnect the electrical cable which powers the motor inside the box (F). By actuating the lever (G, photo 1) remove the motor box from the apparatus. Take off the lower lid (M) of the motor box. Unthread the cartridge filter(O). Insert the new cartridge and close again the motor box by the six special screws.

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The fluorescent tube replacement operations can be carried out by the operator who uses the machine habitually provided that he/she is trained appropriately, after disconnecting the machine power supply.

When the fluorescent tube burns out it is necessary to replace it as described below:


Photo 3-4: Open the cover of the apparatus (A) and block it in the open position by means of the rod (B). Unscrew the handwheel (M), remove the fixing plate (N) and the light fixture (O). Remove the fluorescent tube (P) from its seats and put a new tube in its place by proceeding in the inverse order.

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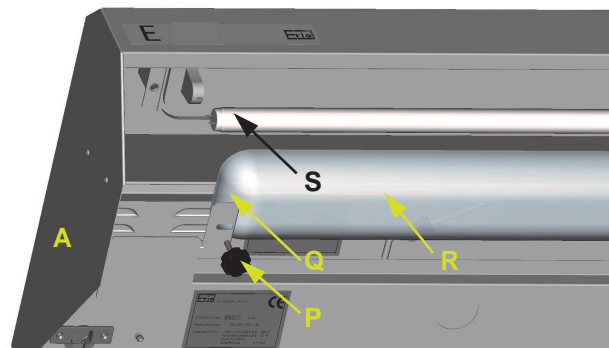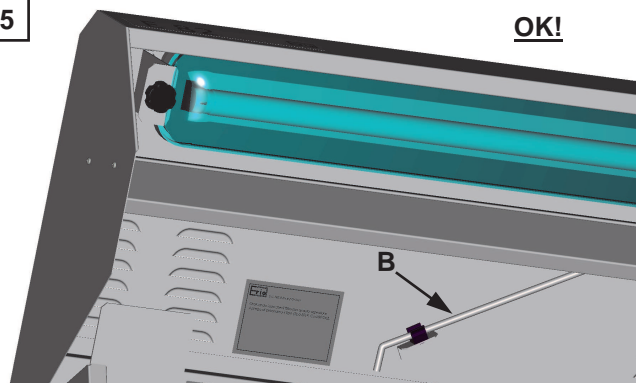
We recommend daily cleaning to remove machining residues and prevent large deposits of material. Use a moist cloth to clean the outside of the machine; do not use alcohol or other solvents.

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2.8 Ordinary maintenance of the machine

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A) Replace the main filter bag when the suction becomes insufficient (2.2).

B) ****The secondary filter cartridge, motor and motor housing should be cleaned every six months, with a small jet of compressed air (2.5).**

C) ****The secondary filter cartridge (2.5), should be changed every 12 months.**

D) **** Only for 110V motors** - After 700 working hours, a check of the motor brushes length has to be carried out. Change the brushes at 900 hours. Following the first replacement, check about every 300 working hours and carry out the replacement if necessary (2.3).

E) Change the fluorescent tube when it burns out (2.6).

**** = The maintenance operations must be carried out exclusively by specifically trained qualified personnel.**

Note: Do not carry out maintenance operations other than those specified in the manual. Any operation not envisaged in the documentation can be hazardous.

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110v 60 Hz

REL0010 800w Motor 110v 60Hz
REL0034 Fuse holder
REL0050 Fuse 5x20 T12A
REL0048 Luminous switch
REL0090 Motor brushes
RFI0042 Main filter bag
RFI0104 Secondary cartridge filter

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The following operations must be carried out by qualified personnel authorized for the purpose, after disconnecting the machine power supply.

For controlling starting and stopping of CS/3000, CS/6000 and AS/TURBO centralized extraction systems, connect the two microswitches (T) located at the sides of the splashguards as indicated in photos 5, 6 and 7.

The control is shown in the OFF position in photo 6 and in the ON position in photo 7.

The connection must be made on the normally open contact of the microswitch.

- The polishing motor supply lead must be connected to the provided connectors inside the connector box.

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Connect the two splash guards (RTB0170 and RTB0172) with Ø50mm (U) flexible pipe to the header RTB0394 (V).

For connections coming out from the wall or the ceiling the drawn air exit pipe (Ø75mm) (Z) can pass through the machine side by opening the appropriate marked hole. If the extraction system is under the floor, the pipe can pass through the machine open bottom and the top, as long as you drill the top with Ø75mm.

Carry out the electrical connection of the control wires, composed of a pair with 2x0.5 mm² section for low voltage, to a CS/3000 or AS/TURBO centralized extraction unit, by means of appropriate connectors/bushes.

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
ATTENTION: The polishing lathe is designed to work at 50% intermittent utilization with a maximum of two hours work combined with two hours of rest. Failure to adhere to this rule may cause the motor to burn out.

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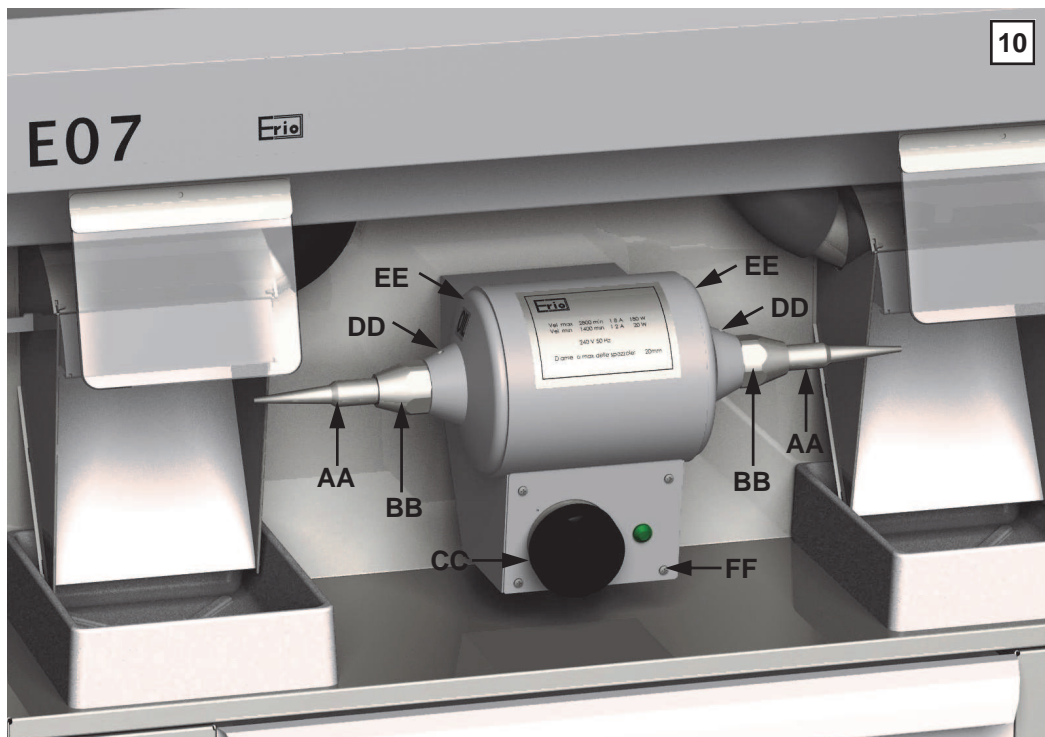Insert the two threaded conical spindles (AA) to the lathe armature having accurately cleaned it of the excess grease. Use particular attention to insert the right threaded spindle (marked with "two rings") on the right armature and the left threaded spindle (marked with "one ring") on the left armature. These spindles are easily removed using the supplied saddle key to turn the extractor nuts (BB). This procedure automatically causes the removal of the spindles. N.B. : It is recommended to perform periodic removal and cleansing of the cones and motor axle to ensure both functionality and durability of the unit.

ATTENTION, before working with this unit use gloves and wear glasses protection.

MOTOR IGNITION: To start the polishing lathe, turn the knob of the ignition switch, located on the front of the base (CC) unclockwise to the position "AVV". This procedure activates the relay switch which starts the unit. The knob automatically will return to position "1"(1400Rpm).

The relay switch, as requested by standards, prohibits the unit from automatically restarting upon the restoration of electrical current after a power failure. In this case, to start the polishing lathe again it has to repeat the above ignition procedure.

VERY IMPORTANT: In the two speed versions, it is recommended to avoid sudden speed increasing and to always start the unit first at low speed (position "1") and then to high speed (position "2"), each time assuring to have turned the ignition knob to the "AVV" position before releasing it. Furthermore, it should always avoid sudden passage from high speed to low speed as this may cause damage to the motor and compromise the lathe's functionality. To avoid it, always turn the machine off from high speed (position "2"+"AVV") to stop (position "0"), wait for ten seconds and then pass to low speed (position "1" + "AVV").

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CLEANING: clean the unit with a damp cloth and never use a steam generator as the vapors may enter the motor and short it.

LUBRIFICATION: The polishing lathe is equipped with special porous bronze bushes instead of ball bearings to ensure durability and silent operation. They require lubrication approximately every 5-6 months. This procedure is easily performed by removing the two screws (DD) and filling the internal reservoir with standard motor oil. Be careful to completely retighten the screws when finished. N.B.: With a new unit or after lubrication, there may be a small loss of oil from the endplates (EE). This is normal; just wipe them clean.

CAPACITOR REPLACEMENT: If the capacitor should burn out, it is easily replaced. First disconnect the electric plug from the wall socket and then remove the four screws (FF) fixing the front panel, taking out the switch and the capacitor. Before disconnecting the fastons and relative electric wires attached to the capacitor, **to avoid electrical shock, take caution to discharge any electrical charges accumulated in the capacitor** by "shorting" the capacitor. Place a screwdriver with an electrically insulated handle assuring the two electrodes at the same time. Replace the capacitor with a new one having the same features like that existent.

FUSE SUBSTITUTION: To change a burnt fuse, it is necessary to remove the fuse cap with a screwdriver, slide out the fuse and replace it with one of equal size (Ø5X20) and power (5A), being carefully to accurately and completely retighten the fuse cap.

IGNITION SWITCH REPLACEMENT: Proceed as in the capacitor replacement by first removing the electrical plug from the wall socket, removing the four screws (FF) fixing the front panel and taking out the switch (CC). Discharging the capacitor as described above. Before disconnecting the switch take note of the exact position/colours of electric wires. Replace it with an identical piece. Position the new ignition switch and close the unit prior to testing.

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The lathe does not start. What to do.

- A.) Control that the electrical plug is plugged into the wall socket.
- B.) Control to have accurately followed the motor ignition instructions
- C.) Perhaps the fuse is burnt. Substitute it as indicated above having first controlled the possible causes as indicated later on.
- D.) Perhaps the capacitor is burned out. Replace it.
- E.) Perhaps the ignition switch is defective. Replace it.
- F.) The relay switch is defective. Call a qualified technician.
- G.) The motor is short circuited. Call a qualified technician.

The fuse is burnt. Substitute it, having controlled and corrected the following possible causes:

- A.) The ignition switch is defective. Substitute it.
- B.) The capacitor is defective. Substitute it.
- C.) The relay switch is defective. Substitute it.
- D.) The signal lamp burnt out. Call a qualified technician.
- E.) The motor is short circuited. Call a qualified technician.

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The E07 machine is supplied in cardboard boxes on pallets. We recommend keeping all the packaging so that you can return the machine to the manufacturer or to the dealer for extraordinary maintenance or repair.

The following items should be found inside the box containing the machine:

- suction unit for E07 machine – E07 instruction manual – polishing motor (if requested) with certificate of compliance.

To limit the possibility of damage and resulting charges for unforeseen repair expenses, we recommend proceeding as follows if the apparatus has to be shipped to the manufacturer for overhaul or repair:

- pack the E07 in its original box. If this is not possible, pack the machine with polystyrene foam protection.

- place the box on a pallet.

- ship it to the manufacturer or local dealer.

The E07 machine must be transported to the point of installation by means of a pallet truck and at least two persons are required for placing it on the worktop.

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6.1 Installation of the protection devices

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To install the two protective screens in Lexan, proceed as described below:

- Having removed the packaging from the machine, open the drawer of the E07 machine and get the two protective screens.
- Raise the cover of the apparatus and block it open by means of the provided rod.
- Fit the two protective screens: the groove in the stainless steel front panel should slot into the front groove in the cover close to the lamp as shown in photos 11-12. Then close the cover.

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. Follow the instructions carefully to ensure correct and safe operation. Keep this sheet for future needs. Contact the supplier for any problems.

. In order to obviate the residual risk of clothes and hair being caught in the brushes the operator shall:

- Operate the machine wearing suitable clothing and avoiding scarves, waist bands and large-sleeved clothing;
- Keep hair tied and /or covered with a hat

. The product may not be altered. Any change nullifies safety approvals and can make the product dangerous.

. Erio S.r.l. shall not be responsible for altered products.

. For safety reasons, those components which have been damaged during operation must be replaced with similar ones before using the appliance again. For maintenance operations disconnect power supply.

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DATA: 10.03.2017 REV.4

PAGINA: 34/39

Erio srl. also declares that the above mentioned machines maintain their compliance with the standard if: they are installed, used and regularly maintained and cleaned according to the recommendations supplied in the instructions for use;

no component has been added, removed, modified or replaced without the express authorisation of Erio srl.;

no other electrical or pneumatic appliance not expressly authorised by Erio srl. has been connected to the machine after delivery.

Responsibility for damage caused by defective products (Dlgs 206/05 and Dir.1999/34/CE) expires 10 years after the date of sale, and then rests with the user.

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DATA: 10.03.2017 REV.4

PAGINA: 36/39

**WARRANTY TERMS ERIO PRODUCTS
General 12-month quality warranty**

All products are guaranteed for 12 months from the date of delivery against material and manufacturing defects. Products found to be defective during the warranty period will be repaired at the Erio works in Zanè (Vi) free of charge or replaced, but no work will be carried out at the Customer's premises.

The analysis of the defectiveness will be carried out by Erio S.r.l. at its works and, if the outcome is positive, Erio S.r.l., at its sole discretion, will replace the defective piece/component free of charge. Shipping and other expenses are not covered by this guarantee in any way.

Products must be returned intact in their original or similar packing and with appropriate outer packaging in order that there will be no damage suffering during transport. Any parts damaged during transport because of poor packaging or other causes not attributable to Erio, and highlighted on arrival of the returned product, shall not be credited in any way.

Please include all your data in the pack along with the reason for the return, and the number and date of the document with which you received the material. The procedure for returning products must be agreed on in advance by fax or e-mail.

Warranty on painted parts:

Erio dental laboratory equipment is built in oven-epoxied galvanised steel plate. The duration of the warranty against corrosion of these materials is 5 (five) years.

The warranty on the painted is forfeited if pickling acids, strong bases or acids and/or corrosive materials are used on them.

Warranty on stainless steel parts:

The bench tops in AISI 304 stainless steel are built with certified materials and the guarantee against corrosion is for 5 years (five).

The warranty is forfeited when:

- A) use is made on the stainless steel parts of materials not permitted according to the merchandise tables attached to the tops.
- B) there are phenomena of electrolysis and electro-corrosion caused by special equipment or failure to earth the stainless steel tops.
- C) episodes of "pollution" occur with iron particles (filings) deposited on the stainless steel surfaces.

Warranty on the electric and electronic parts:

The electric and electronic equipment is guaranteed against manufacturing defects for 1 (one) year. The warranty expires when repairs or tampering are carried out by personnel not authorised by Erio srl in Italy, or by an Official Erio Dealer in other territories.

The warranty for electric motors is valid for up to 2500 hours of effective use or for one year, whichever is the more severe.

The warranty for 110V electric motors is valid for up to 700 hours of effective use or for one year, whichever is the more severe.

Starting date:

The warranty on Erio s.r.l. products runs from the date of installation of the material in the user-customer's laboratory. This date must be witnessed by the waybill which accompanied the goods or by the date of the test report issued by the personnel responsible for assembly.

Exclusions:

Damage caused by improper use of the products is not covered by the warranty, and no responsibility is accepted for damage due to such improper use.

It remains understood that the guarantee is forfeited at the moment products supplied by Erio S.r.l. are used outside the field of application for which they were designed, that is, outside dental laboratories.

The warranty on electric motors of suction groups expire if you use non-Erio filters.

Accessories are not covered by the warranty if they are not installed on Erio units for which they were expressly designed.

The guarantee is subordinate to correct installation which must be carried out by personnel trained and authorised by Erio S.r.l. or by the official Erio Dealer for the territory in question.

Such personnel is trained through support or by means of training courses held in the Erio works.

We decline all responsibility for consequences or damage deriving from the goods supplied by us.

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